**SPORT-ELEC®**

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## I. INTRODUCTION

You have just purchased a product designed and manufactured by SPORT-ELEC®. We thank you and invite you to read the user manual before starting your session.

SPORT-ELEC® is a french designer, manufacturer and distributor of consumer and professional electro-stimulation devices for over 20 years. Based on the latest innovations from the Sport-Elec research office, COMBI optimizes the efficiency of electro-stimulation and sports sessions.

The COMBI is a set of 6 accessories including:

- 1 vest + 8 electrodes
- 1 belt + 4 electrodes
- 2 thighs + 4 electrodes
- 2 armbands usable on the calves + 4 electrodes

## II. COMBI ACCESSORIES

The 6 Combi accessories allows you to work out several muscle groups at the same time and very precisely with any electro-stimulation device.

Each accessory comes with adhesive electrodes that adhere to the skin and have to be scratches inside the accessories. Position the electrodes where you want to work out according to your morphology and body size. Connect the electrodes onto the accessories and to the electro-stimulation device. Moisten the electrodes with water and adjust the accessories around your body.



- ▶ Usable with any electro-stimulation device \*
- ▶ Fits men and women size S to XL
- ▶ Position electrodes wherever you want
- ▶ Accessories can be used for a session of electro-stimulation in static position or during a sport training

\* Adapter available for sale on the website [www.sport-elec.com](http://www.sport-elec.com).

The COMBI is adaptable with any electro-stimulation devices using plugs of 2 mm diameter.

## III. CONTENT OF THE PACK

An electro-stimulation vest with 8 adhesive electrodes to position inside the accessory and stimulate the abdominals, pectorals, dorsal, lumbar or trapezius muscles

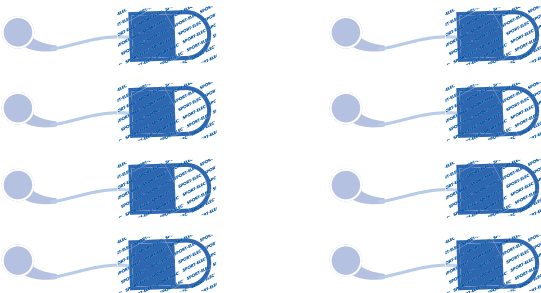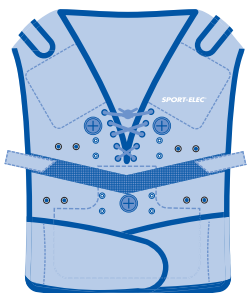
## IV. USE

### 1 Position the electrodes inside the accessory


Position the electrodes inside the accessory following the examples below to stimulate the front body muscles - Abs, Pectoral - or the back body muscles - dorsal, lumbar, trapezius.

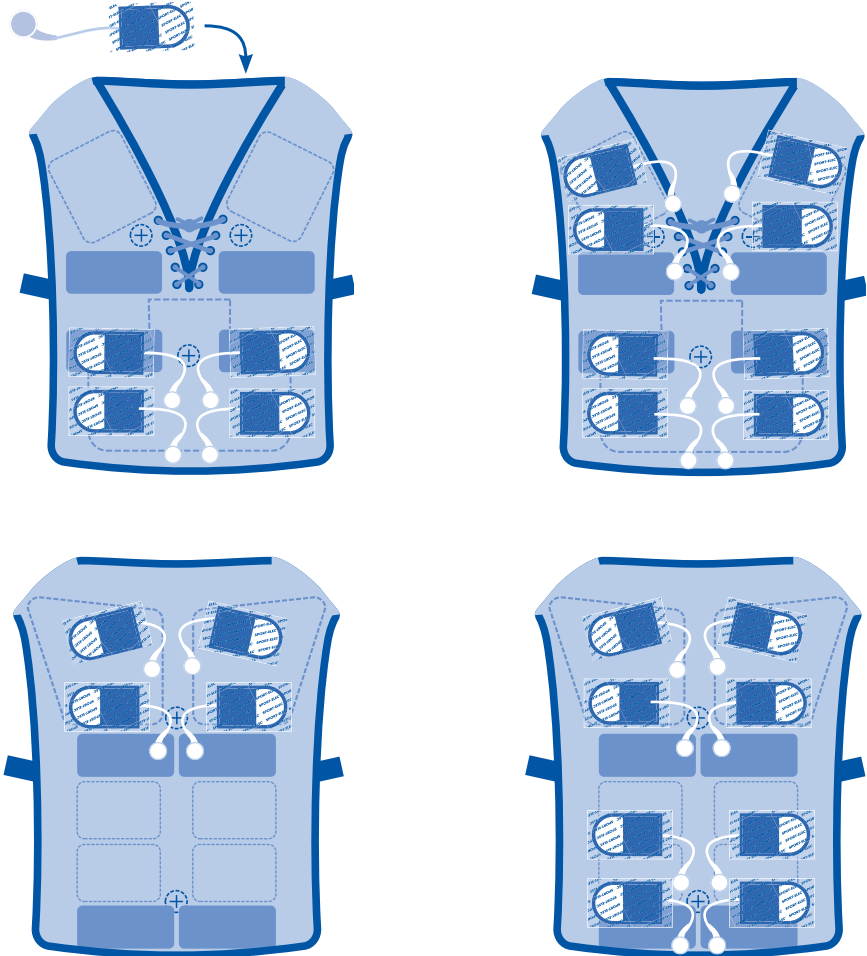
## 2 Pass the electrode wires out of the accessory

Pass the wires of the electrodes in the slot of the accessory

1

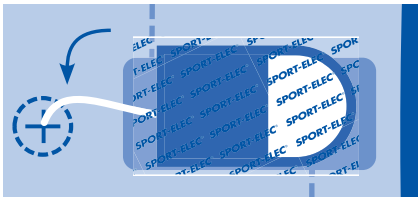2

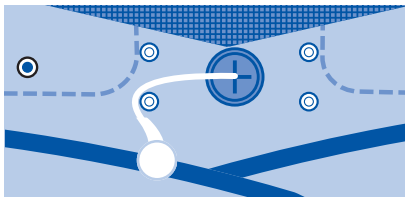
## 3 Clip the electrodes to the snaps of the accessory

Connect each electrode to the outside of the accessory

1


2

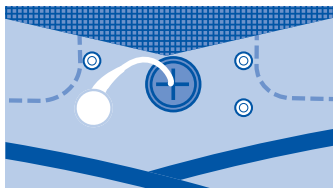3


4

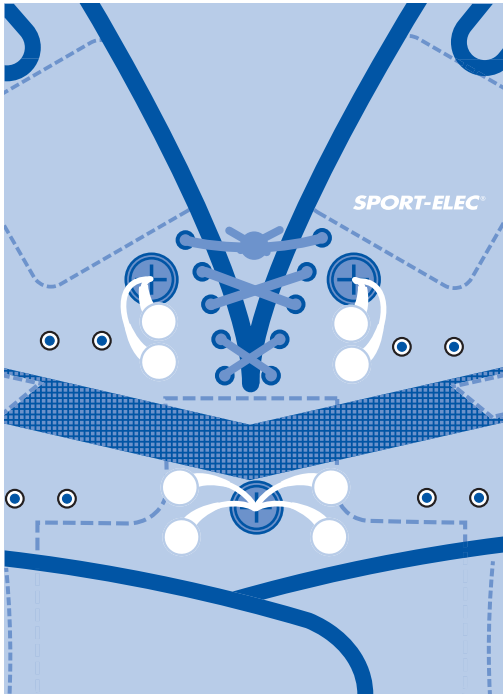
#### 4 Connect the generators of your electro-stimulation device onto the accessory

Clip the 2 generators on the external side of the accessory. For example with MultiSport Pro or FreeAction devices.
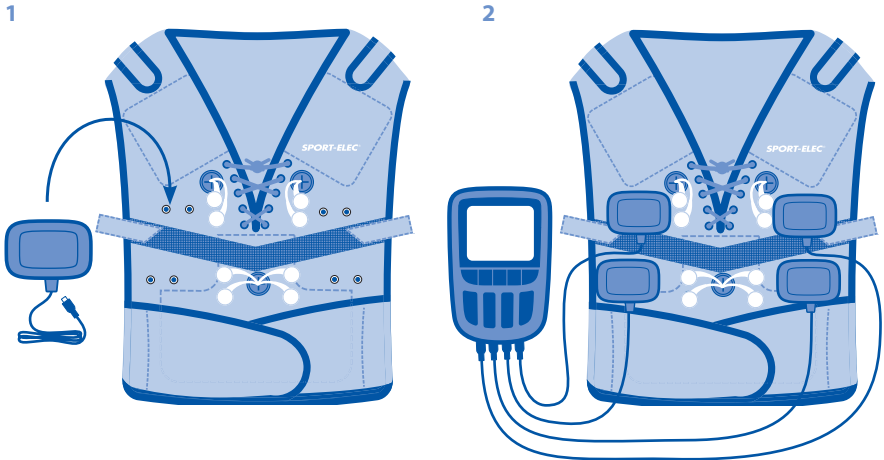
#### 5 Remove the plastic from the electrodes and slightly moisten with water

Moisten electrodes to improve current conduction and comfort.
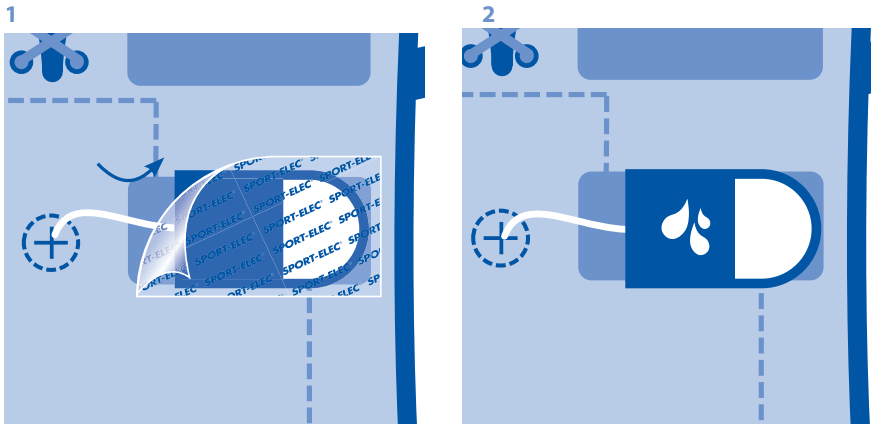
## 6 Put on the vest

Put on the vest to stimulate the muscles of the front of the body (Abs, Pectoral) or the back of the body (dorsal, lumbar, trapezius).

## 7 Start your electro-stimulation device, choose the program and increase the power.

Enjoy !

# V. EXAMPLES OF COMBI USE WITH SEVERAL ACCESSORIES

## 1 Electro-stimulation device using 2 generators

If you have an electro-stimulation device using 2 generators you can work out 2 muscle groups simultaneously:

- ▶ Arms (biceps or triceps) using armbands the calves (internal and external twins) using the armbands
- ▶ the abs using the vest
- ▶ the pectorals using the vest
- ▶ the dorsals using the vest
- ▶ the shoulders using the vest
- ▶ the gluteal using the belt
- ▶ the thighs using the thighs accessory

## 2 Electro-stimulation device using 4 generators

If you have an electro-stimulation device using 4 generators you can work out 4 muscle groups simultaneously:

- ▶ the pectoral and the abs with the vest
- ▶ the back and shoulders with the vest
- ▶ arms and thighs using armbands and thighs
- ▶ the inner and outer sides of the thighs with thighs accessory
- ▶ the abs and gluteal with the vest and the belt
- ▶ the abs and thighs with the vest and thighs
- ▶ gluteal and thighs with the vest and thighs

# VI. TECHNICAL CHARACTERISTICS

## 1 The accessories

Material: black polyester jersey fabric with velcro straps.

Connectors: snaps.

## 2 The electrodes

Material: conductive hydrogel with velcro strips on the outside.

Connectors: wired with snap.

# VII. CARE AND CLEANING TIPS

## 1 The accessories

- ▶ Once your session is over and the electro-stimulation device disconnected from the accessory, use a damp sponge or a dry cloth to gently clean the fabric if needed.
- ▶ Do not machine wash • Do not dry clean • Do not immerse in water • Do not use bleach.

## 2 The electrodes

- ▶ At the end of the session we recommend moistening the electrodes again before putting the plastic back on.
- ▶ To extend the adhesive electrodes usage, remove it from the accessory and store it in a dry place. Some users recommend that you keep the adhesive electrodes in the crisper drawer of your refrigerator.

## VIII. DO AND DON'T

- ▶ Follow the precautions and contraindications of the electro-stimulation device that you connect to the COMBI.
- ▶ DO NOT USE ON THE HEART AREA OR ON THE GENITALS, DO NOT USE IF YOU ARE PREGNANT, DO NOT USE IF YOU HAVE HEART CONDITION OR WEAR A PACE MAKER.

## IX. SECURITY CONTROLE

Every 24 months, check the following:

### 1 Check if there is a superficial degradation of the material

### 2 Check the reliability of the accessories

No tears, no defects in materials.

## X. WARRANTY

The accessories are for personal use only. Sport-Elec does not take it back or replace it.

## XI. RECOMMENDATIONS

When the device is out of order and to be destroyed please be sure to follow national regulations. Please use the convenient facility or location to through the device away in order to ensure its safe disposal in an environmentally friendly manner.

## XII. TRANSPORT AND STORAGE

- ▶ Storage temperature between + 10 ° C and + 40 ° C
- ▶ Relative humidity between 30% and 75%
- ▶ Pressure 700 hPa.

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